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Thank you for purchasing our product "Exidéal libra".  
Please read this Manual carefully before using this product.  
Please check the QR code on the right for the full text of the Manual.

### Before Using

- Please connect the included USB charging cable to your USB power source to charge the battery.
- \*5V1A USB power source recommended.
- When irradiating this product to the face, be sure to use the dedicated goggles.
- If the goggles slip or slide, attach the attached nose pad to the nosepiece and adjust the goggles.

### How to attach the goggles

Insert the handle of the dedicated goggles into the socket on the inside of the controller by sliding it from the tip to the back.  
When installing or removing goggles, be sure to hold the base of the goggle slot firmly with your fingers as shown in the illust ration.

- \* Forcing it in or out may cause malfunction or damage to the goggles.
- \* When not using, be sure to remove the dedicated goggles and store it.


### Program List

You can select of programs according to your skin condition and troubles. For continuous use, 12 minutes or less is recommended.

LED Program	Features	Time	Area of Use
Texture/ Resilience	For age-appropriate care, aim for textured skin with tense.	4 Min	Face Neck
Adolescent skin	Support skin that becomes unstable due to disturbed life rhythms and dietary habits to smooth skin.	4 Min	Face Neck
Scalp care	Age-appropriate care to keep the scalp and hair in good condition. Used to care repeated coloring and UV-exposed scalp.	4 Min	Scalp

### How to Use

- Long press the Power/Program Button to turn on the power.  
→ The irradiation of [Texture/Resilience] will start automatically with beeping sound.
- Each time you press the Power/Program Button, the program will change in the following order:  
[Adolescent skin] → [Scalp care] → [Texture/Resilience] → [Adolescent skin]
- To power off the program while it is operating, long press the Power/Program Button again.  
→ The power will turn off with beeping sound.
- \* The program cannot be paused.
- \* Switch the button during irradiation, the new program will be started for 4 minutes.

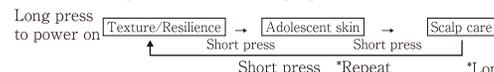- When the 4 minute program ends, the power turns off automatically with beeping sound.
- When the battery is runs out, the indicator lamp will flash 5 times with beeping sound. Charge the battery as soon as possible.  
\* While charging, the main body cannot be operated.

※ Please keep this warranty in an accessible place for easy reference.  
※ Misuse and incorrect storage of this product are not covered by the warranty, and repair is subject to charge even during the warranty period.