

Warranty certificate Exidéal ovo

This warranty guarantees that if a failure occurs within the following period from the date of purchase, free repair will be served as described below.

The purpose of use of the personal information you have written is described below. Please refer.

| | | |
|----------------------|-------------------------|--|
| Stock Number | EX-HA01PG | |
| Serial Number | | |
| Warranty Term | One year after shipment | |
| Purchase Date | | |
| Customer information | ZIP | |
| | Address | |
| | Name | |
| | TEL | |
| Distributor | | |

【Regulation for free repair】

1) If the product fails during the warranty period under the correct usage according to the instruction manual and the label attached to the main body, free repair will be served based on the contents of this warranty card.

● Please bring the warranty when you need the free repair.

Shipping charges for all warranty repairs are at the owner's expense when you send it to our customer service or original distributor.

● Please contact our customer service if you do not ask for free repair at the original distributor.

2) Please contact our customer service about free repair when you move to another place.

3) Please contact our customer service when you can not receive free repair service at original distributor due to the product being given as a gift.

4) This warranty does not apply:

● The damage caused by misuse, repair, alteration,

● The damage caused by relocation of storage, transport, drop after purchasing

● The damage caused by liquid contact, fire, earthquake or other external cause, to incorrect use with abnormal voltage or an undesignated power source (voltage, frequency)

● Shipping damage or damage caused by mishandling by the freight carrier

● The damage caused by use with other than home use (Example: business purpose)

● Not provide this warranty

● Not provide proof-purchase day, distributor name, your name and reformatted in this warranty

5) Consumable items are not covered by the warranty.

6) This warranty is only valid in Japan.

7) We do not reissue this warranty.

■ Purpose for usage of the personal information

※ We will maintain and use customer information (copy of warranty) for free repair support and safety inspection within the warranty period.

※ This warranty guarantees that if a failure occurs within the following period from the date of purchase, free repair will be served. It does not limit the customer's legal rights and obligations to the person responsible for warranty and other business partners-, issuing this warranty. Please contact our customer service when you have questions about repair after the end of the warranty period.

※ For after-sales service after the warranty period has expired, refer to " After purchase servicing " in the manual.

Haslux Co., Ltd.

ZIP: 106-0032

29F ARK HILLS SENGOKUYAMA MORI TOWER, 1-9-10 Roppongi, Minato-ku, Tokyo, Japan

TEL: +81 03-6206-6188 FAX: +81 03-6206-6289


<https://www.exideal.jp>

Warranty attached

Manual / Warranty



| | |
|---------|---|
| P.1-P.2 | Safety Notes |
| P.3 | Features of Exidéal ovo /Before using |
| P.4 | Items included and name of each part / Preparing |
| P.5 | Name of button and function / Explanation of display / How to operate |
| P.6 | Mode List / How to install and remove the dedicated cotton |
| P.7 | Tips for beauty care / After using |
| P.8 | Troubleshooting |
| P.9 | After-sales service / Q&A |
| P.10 | After purchase Service / Specifications |

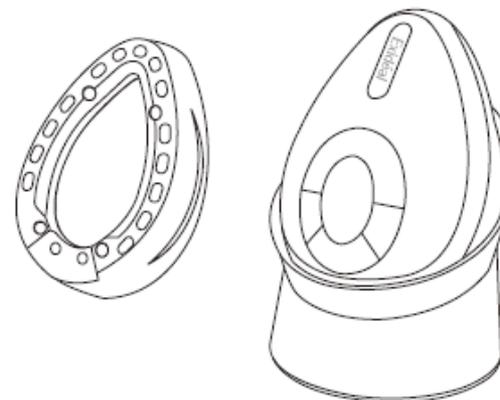

Thank you for purchasing our product "Exidéal ovo".
Please read this manual carefully and use our product.
The warranty is included on the back side of this manual. Please keep it safe.

Haslux Co., Ltd.

Please read this instruction manual carefully and use the product safely and correctly.

- These notes provide instruction on the safe and correct use of our product and to prevent injury and damage.
- Below are the classifications for these notes;

Caution classification

⚠ Danger This symbol indicates an imminently hazardous situation which, if not avoided, may result in death or serious injury.

⚠ Warning This symbol indicates that the user risks of death or serious injury if not avoided.

⚠ Caution This symbol indicates a potentially hazardous situation, which, if not avoided, may result in personal injury or damage to the equipment.

Explanation of Figure

- ⊘ This symbol is displayed for acts prohibited.
- This symbol Indicates a mandatory action that must be followed.

⚠ Danger

Persons who use a heart pacemaker or electro-medical apparatus that may cause malpractice. Please do not use this product.

⚠ Warning

When using this product, do not look directly into LED light.

Do not use this product for any other purposes than making your face and body beautiful, or do not use the broken product. If you have any concerns about your skin and physical condition when using, please consult a doctor before using.

For those who have any skin problems listed below, please do not use this product. It may worsen skin problems.

- Photosensitivity
- Those who have skin disease or medical history of skin problems
- Dermatitis (Atopic dermatitis, Contact dermatitis, excessive sunburn, etc.)
- Metal allergy
- Areas with purulence due to acne and pimple
- Eczema, swelling, and warts (including virus)
- Purpura
- Melasma

For those who have any disease listed below, please do not use this product. It may cause unexpected accidents or skin and body problems.

- Acute disease
- Blood disease (hemorrhagic disease, hemophilia, etc.)
- Those who are being under treatment or being followed up for abnormal blood pressure · Febrile disease · Malignant tumor
- Epidemic (Infectious diseases) · Visceral disease
- Facial neuropathy
- Head and ophthalmologic sequelae due to traffic accidents
- Those who have any abnormality in bone (osteoporosis, fracture, etc.)
- Acute pain disease (sprain, muscle strain, etc.)

⊘ Do not use this product on the skin and the areas that follow below. It may cause unexpected accidents or skin and body problems.

- Sensitive skin
- Skin Paresthesia (especially those who cannot recognize heat)
- Bulging mole
- Cut
- Surgical site due to accidents and orthopedic surgery
- On implant (implant surgery metal / silicon / plastic, etc.)
- Around the eyelids, eyeballs, and hollows around eyes
- Adam's apple
- Mouth and genital area, other mucous membrane
- Thin skin area (immediately above the bone)
- Around the heart

⚠ For those who have any of the following , please do not use this product. It may cause unexpected accidents or skin and body problems.

- Those who have epileptic seizures that respond to blinking lights
- Allergic constitution
- Those who have telangiectasis due to long-term use of steroid hormones or due to liver dysfunction
- Those who are during menstruation
- Those who are pregnant, one month after childbirth, and during nursing
- Those who are under dental treatment, or temporomandibular disorders
- Those who are in a growth process
- Those who are after drinking alcohol or taking medicine (especially medicine that induces sleep)
- Those who are very tired, who do not feel well, and who need rest
- Those who are with body temperature of 38°C or higher (fever period)

Do not use this if there is any damage and abnormal heating in the cord and the dedicated AC adapter, or if it is loosely inserted into the outlet. There is a risk of fire, electric shock, short circuit, and burns.

Do not hurt, damage, modify, forcibly bend, pull, twist, place a heavy object on the cord. There is a risk of fire, electric shock, short circuit, and burns.

Be sure to use with AC 100-240V outlet. Also be sure to use the dedicated AC adapter for this product. If not, it may cause fire, electric shock, and failure.

Remove precious metals such as earrings or necklaces for use. It may cause accidents or problems.

⚠ Warning



Do not store this product in a bathroom or in a humid place, and do not place it on where is subject to water (such as a washbasin). There is a risk of fire or electric shock due to deterioration.

Do not throw this product into fire or heat it, charge and use it , or leave it alone at a high temperature. If not, it may cause heating and fire, burst, because this product has rechargeable battery built-in.

Do not use this product near any flammable object (alcohol, thinner, benzine, liquid remover, spray, etc.) If not, it may cause fire or electric shock, or injury.



In the following cases, please stop using this product, turn the power off, and unplug the dedicated AC adapter from the outlet. There is a risk of fire or burns due to abnormal output.

- Smoke from this product
- Abnormal heating from this product
- Abnormal noise from this product
- Odor from this product



Never modify, disassemble or repair this product except by repair service engineers designated by us. There is a risk of fire or injury due to abnormal operation.



Do not soak this product in water. There is a risk of electric shock or short circuit.

Be sure to turn the power off when not in use or when cleaning this product. It may cause electric shock, short circuit, or fire.

Unplug the dedicated AC adapter from the outlet except when charging. It may cause electric shock or earth leakage fire due to insulation deterioration.

Do not connect or disconnect the dedicated AC adapter with wet hands. Also, make sure there is no water on the charging terminal before charging. There is a risk of electric shock or short circuit.



This product contains a very strong infrared irradiation pulse program, infrared remote control devices such as television or air conditioner may not work. In the event that malfunction or any situation where remote control cannot be used arises, please turn off the power of this product.

Keep this product out of reach of children. If not, it may cause accidents and problems.

Do not use with pets or young children, and do not let them play with this product. If not, it may cause accidents or injury.

For those who use this product for the first time or who have sensitive skin, shorten usage time to use it while checking the skin condition. If not, it may cause accidents and problems.

Follow the instruction manual for usage time and frequency. If not, it may cause accidents or problems.

If you feel pain or discomfort, stop using this product. If not, it may cause accidents and problems.

Do not keep applying this product on the same area and be sure to move it on the skin. If not, it may cause burns or low temperature burns.

⚠ Caution



On the day of hair removal or shaving process, do not apply this product on the same area. There is a risk of swelling and redness.

Do not drop, throw, and give a strong shock on this product. If not, it may cause accidents and problems.

Do not place anything on this product.

When unplugging the dedicated AC adapter, be sure to hold the adapter plug and not to pull on the cord.



Not to store this product in the following places.

- Place exposed to direct sunlight
- Place with high humidity
- Place with a lot of dust
- Place where there is a risk of ignition
- Place exposed to wind and rain
- Place at high temperature or below zero
- Place where chemicals or gas are generated
- Place with violent vibration

Do not wipe this product with chemicals such as alcohol, thinner, benzene, alkaline detergent, bleach, neutral detergent, or disinfectant. If not, it may cause malfunction such as deformation, discoloration, crack, or damage.

※Please keep this warranty in an accessible place for easy reference.

※Misuse and incorrect storage of this product are not covered by the warranty, and repair is subject to charge even during the warranty period.

Features of Exidéal ovo

For the best skin ever, with Double-power of Exidéal's original LED program × 4 esthetic functions.

Multi-functional handy model made by Exidéal, a pioneer of LED facial beauty device. It is 5in1 beauty device, with a unique combination of specific lights (wavelength) that combine LED programs allowing you to give treatment according to the purpose, and 4 esthetic functions leading to the selected beauty skin. 8 patterns* of care for each problem can support your skin care throughout a year to fulfill your ideals.

① Double-power of LED program × esthetic function

We can realize efficient maintenance with Double-power by combining LED program and esthetic function that suits your skin problems, because of the amount of years we have researched LED.

② Small LED with high luminance and high output

By adopting chip type LED and calculated layout design, we have achieved to equip with high luminance and high output LED (55pcs of 4 wavelengths) even with compact handy facial beauty device. You can get real beauty care at home that you've been able to get only at beauty salon.


③ Circle optical design

A unique optical design that increases the irradiation efficiency by 120 to 130%, with the calculated LED arrangement and angle, and the mutual reflection of the mirror-finished reflector (mirror panel) and skin, enable more efficient photo rejuvenation.

* Compares with a beauty device that is not equipped with mirror reflector

④ Body care also possible


You can use not only for your face, but also for all parts of your body. Mode 01(EMS) is recommended for upper arms and thighs which are especially concerned.

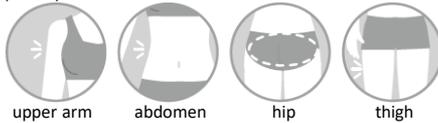


*Including esthetic function ON/OFF

Items included and name of each part


● Main body

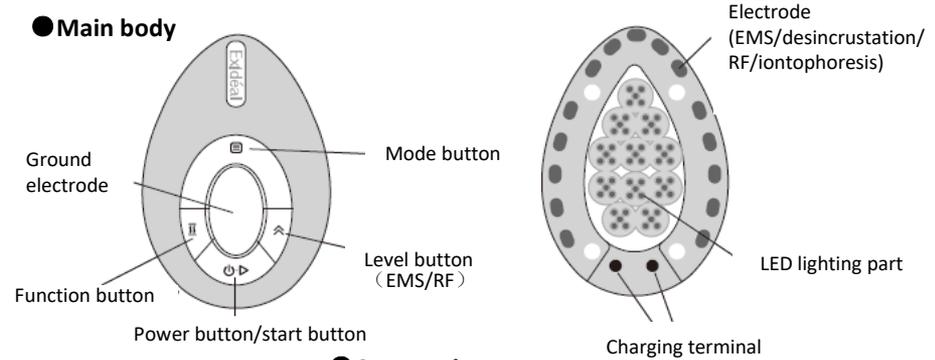

● Charging stand

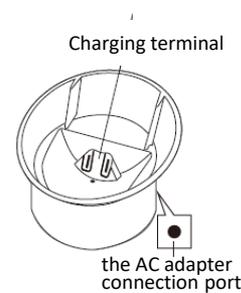

● Accessories

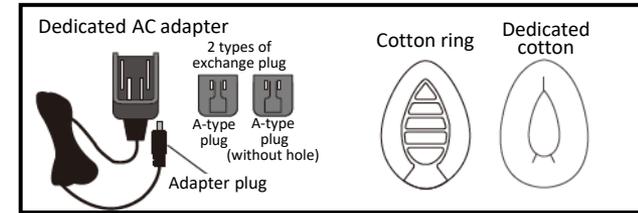

This product comes with two types of plugs that can be used in each country. If it does not match your outlet, please use a conversion plug, which matches your outlet.

*For use in Japan : A-type plug
 For use in China : A-type (without hole)

Before using

◇ Frequency of use

Use within 24 minutes in total a day.
 Do not use esthetic function EMS/RF more than once a day.

◇ About skin

Before you use, be sure to perform a patch test (skin test) for each cosmetic you use. (P3)
 *This is a simple test to check any skin problems in advance.
 Use this product with your bare skin (after removing make-up and facial washing).
 Use this product while checking your skin condition.
 *If you have any skin problem, please stop using this cosmetics with the product and consult a doctor.

◇ Accessories

Before using, please remove the accessories around your face and neck such as earrings.

Patch test (skin test)

When you use this product and cosmetics for the first time, be sure to perform a test to prevent skin problems in advance.

- 1 Wash the area (inside the arm) to be performed a patch test with soap, and cleanly wipe with towels.
- 2 Apply cosmetics (lotion, beauty essence)at inside part of the arm, and select the desincrustation (mode 02) or iontophoresis (mode 4).
- 3 Apply the electrode of the product to the same area above, keep moving for about 3 minutes .

If there is no abnormalities within 48 hours, you can use this product.

- * Be careful with any abnormalities on your skin when you use.
- * Do not use the product if it does not suit your skin condition.

Preparing

Charge

When you use this product for the first time or when the battery is exhausted, charge the battery on the charging stand.

- 1 Insert the device plug into the AC adapter connection port on the charging stand.
- 2 Connect the dedicated AC adapter to the outlet.
- 3 Place the main body on the charging stand. —blinking (green) · sound (beep twice)
- 4 Charging is completed in about 120 minutes. — lighting on (green) · sound (beep three times)

○ Approximate charging time : about 120 minutes

○ Startup time : about 24 minutes

* Differs depending on the mode used and operating environment.


- While charging, the main body cannot be operated.
- Even if it is charged more than 120 minutes, the main body is not hindered.
- Time to charge

* If the battery is exhausted while in use, the power button/start button will blink rapidly (yellow) and beep sound is emitted for about 5 seconds. When the battery runs out, the power will be turned off with beeping sound.

<Warning > · Be sure to use the dedicated AC adapter. Or it may cause accidents or failure.

- Connect the charging stand and adapter plug first before placing the main body on the charging stand, so as not to occur short circuit.
- Before charging, please wipe off the lotion and beauty essence on the main body.

Name of button and function

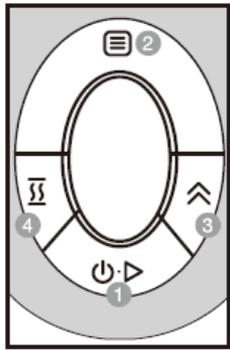

- 1 Power button/Start button**
Long press the button to turn the power on.
Press the button briefly while the power is on will start each mode.
When you want to stop using, long press the button to turn off the power.
- 2 Mode button**
Press the button can select the mode (01→02→03→04→01).
When you want to switch the mode while in use, select the mode again and press the power button/start button to start new mode.
- 3 Level button (EMS/RF)**
When you press the button when using mode 01 and 03, you can select strength of EMS/RF.
You can also switch the strength during use.
- 4 Function button**
When you press the button, you can select ON/OFF of esthetic function (EMS/desincrustation/RF/iontophoresis) in each mode. You can switch ON/OFF during use.

Explanation of display

| | | | | | | | | | |
|----------------------------|--|----------|------------|---------|---------|----------|------------|---------|---------|
| Power button /Start button | Light off: Power off Light on in yellow :Power ON Blink in yellow : Low battery Light on in green : Charge completed Blink in green : Charging | | | | | | | | |
| | | Mode 01 | | Mode 02 | | Mode 03 | | Mode 04 | |
| Mode | Light on: Get ready Blinking: Operating | | | | | | | | |
| Level (EMS/RF) | Light off: Weak Light on: Strong | Weak | Strong | — | | Weak | Strong | — | |
| Function | Light off: OFF Light on: ON | ON | OFF | ON | OFF | ON | OFF | ON | OFF |

How to operate

Common operations of all modes

- Remove the main body from the charging stand.
- Long press the power button/start button will turn the power on with beeping sound.
- Select the mode with mode button.
- Select ON/OFF of esthetic function (EMS/ desincrustation /RF/iontophoresis) of each mode with function button.
- Select the strength of EMS/RF with level button. (Only for mode 01 and 03)
- Press the power button/start button to start operation.
 - When using desincrustation (Mode02) or iontophoresis (Mode 04), please hold the main body in the way that your hand can touch the earth electrode.
 - When you switch the mode while in use, select the mode again and press the start button briefly. Then new mode will start.
- Each mode is completed within 6 minutes* from the start, and the power turns off with beeping sound.
 - If you want to stop using the device in less than 6 minutes*, press and hold the power button/start button to turn off the power with a beeping sound.
- If the power button/start button blinks during use and you hear a beeping sound for about 5 seconds, the battery level is low. Please turn off the power and charge the battery.

*The end time may vary depending on the operating environment

Mode List

You can select from 4 modes for your skin condition and concerns. Also, you can turn off the esthetic function to suit your mood, so you can care with only LED program.




Cotton mark is to recommend of using the dedicated cotton.

- Blue : Adolescent skin / T-zone, etc.
 - Yellow : Texture / resilience, prevention of rough skin
 - Red : Aging skin* / dry skin *care according to age
 - Near infrared ray : resilience / skin environment
- *Each mode consists of unique program that combines multiple wavelengths.

| Mode | Time | Features LED program x esthetic function | Skin trouble · concerns |
|-----------------------------|--------|---|--|
| 01 Texture / resilience | 6min.* | LED program (texture/resilience) × EMS Your muscles are stimulated and strengthened by EMS function, and LED program makes your skin resilience. | Lack of skin resilience, disorder of texture, care according to the age |
| 02 Adolescent skin | 6min.* | LED program (adolescent skin) x desincrustation With desincrustation function, it absorbs dirt that cannot be removed by washing, and LED program leads unraveling smooth skin. | Care the skin that is affected by disturbed lifestyle rhythms and dietary habits. |
| 03 Rough skin / dry skin | 6min.* | LED program (rough skin / dry skin) × RF RF function warms and conditions the skin with frictional heat. LED program leads to moist and gentle skin. | Rough skin /dry skin concerns. Irregular lifestyles and seasonal changes fluctuate the skin. |
| 04 Bright skin | 6min.* | LED program (bright skin) x iontophoresis Iontophoresis function promotes penetration of beauty ingredients into the stratum corneum of skin, and provides it moisture. LED program leads to bright skin like pottery. | Transparent** skin and dullness concerns. Need of aftercare of UV ray. |

*The end time may vary depends on the operating environment.
**Due to skin waste.

How to install and remove the dedicated cotton

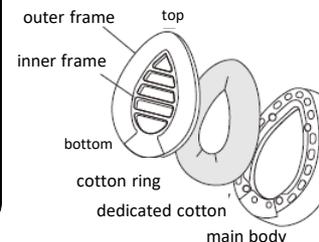

● How to install


- Put a sheet of dedicated cotton on the electrode part.
- Push the inner frame of cotton ring from bottom to top.
- Push the outer frame of cotton ring from bottom to top.

● How to remove

Remove from under (logo part) part of the cotton ring.

<Caution>

If you hold the top of outer frame to remove the cotton ring, it may break.

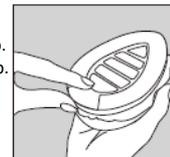


Tips for beauty care

How to hold the body

Hold the body as image on the right.

When using desincrustation(mode 02) or iontophoresis (mode 04), hold it the way that you can touch the ground electrode.

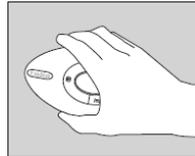

Mode01 Texture/resilience

When you use EMS, make sure that the entire electrode part touches your skin, and slowly moving the body while using. (Reference : 2cm per second)
EMS gives direct electrical stimulation to the muscles, and it is characterized by sharp stimulus. For those who use for the first time, please use it in "weak" mode. There are individual differences in how to feel EMS, so please select the level which you feel comfortable.
Do not use around eyelids or eyes.

Mode02 Adolescent skin

Mode 04 Shiny skin

When using desincrustation(mode 02) or iontophoresis (mode 04), use the cotton ring to attach the dedicated cotton (P6). When using it with esthetic function OFF, the cotton is not necessary.
After putting the cotton on, apply your own lotion or beauty essence to penetrate through entire cotton.

◆Recommended cosmetics

- desincrustation : cleansing lotion, lotion
- iontophoresis: contains Vitamin C, Vitamin C derivatives

After use, remove the cotton ring and dispose the used cotton (P6), and wipe the LED lighting part and the electrode part of main body with a clean cloth. Wash the cotton ring with water and dry it well.

Mode03 Rough skin/dry skin

When using RF, use it that the electrode part touches the skin while moving the main body slowly. (Reference : 2cm per second)

Keep applying it to one area may cause burns or low temperature burns.

If you feel hot, lower the RF level and move it a little faster.

Do not use around the eyelids or eyes.

After using

- Wipe the dirt on LED lighting part · electrode part with a cloth which squeezed tightly after dipping in water.
 - Do maintenance with the power off.
 - Do not use the chemicals such as alcohol, thinner, benzine, alkaline detergent, bleach, neutral detergent, sterilizer, or disinfectant.
 - Do not leave this product with lotion and beauty essence remained.
- Wipe off the dirt on the main body with a dry cloth.
 - This product is not waterproof, so do not wash it with water.

◆How to store

- When you do not use this product for a long time, please clean the main body and store it.
- When not using, please turn off the power of the main body, and pull the dedicated AC adapter from the charging stand and the outlet before storing.
- Do not store this product near a stove or heater, in hot and humid place such as bathroom, place exposed to direct sunlight, or in unstable place.

Troubleshooting

| Symptom of trouble | Possible causes | How to deal with |
|---|--|---|
| Cannot Power On | Power button · start button is not pressed. | Long press the power button · start button (about 3 seconds) to turn the power on. |
| | Insufficient charge. | Please charge the battery. (P.4) |
| The power turns off during use. /Can only be used for a few minutes even after charging. | Insufficient charge. | Please charge the battery. (P.4) |
| | Auto-off function worked. * If no operation for 6 minutes, the power will be automatically off. | Press the power button /start button again to turn the power on. |
| | The mode in use has ended. | When each mode ends, the power turns off automatically. |
| | Battery life is up. | Battery needs to be replaced. Please contact our customer service. |
| Charging does not start. | Dedicated AC adapter is not inserted into the outlet. | Insert the dedicated AC adapter plug into the outlet. (P.4) |
| LED does not light on. /Partly does not light on. | Internal failure | Restart the machine. If the problem is not resolved, please stop using and contact our customer service. |
| | The near-infrared light is colorless and transparent, and it is invisible even during irradiation. | |
| Desincrustation/iontophoresis is weak, cannot be felt. | | Ion is almost not felt by the body. Also, there are individual differences in how to feel. |
| The mode does not end even after 6 minutes. | The end time may vary depending on the usage environment. | |
| | | Long press the power/start button (about 3 seconds) to turn the power off. If the problem is not resolved, please stop using and contact our customer service. |
| EMS electrical stimulation is weak, no stimulus is felt, the warmth of RF is not felt. | Insufficient charge. | Please charge the battery. |
| | Multiple electrodes do not touch the skin. | Let it in close contact to the skin. |
| | Lack of skin moisturizing. | If the surface of the skin is dry, it may be difficult to feel EMS/ RF. |
| The main body is heating. | Continuously used for long time | Continuous use for long time may make the main body slightly hot. Turn the power off and wait for some time before using. If the problem is not resolved, please stop using and contact our customer service. |
| During charging, blinking green light and sound do not stop. | The temperature sensor is broken or disconnected. | Remain the main body from the charging stand. Please contact our customer service since it may be broken. |

Caution

Please contact our customer service if this troubleshooting guide does not apply your problems.

After-sales service

◆ Rechargeable battery

- Rechargeable battery is a consumable item. (Not covered by warranty)
- If the operating time is extremely short after charging, it is considered that the rechargeable battery has reached the end of life.
- If you want to replace the rechargeable battery, please contact the place of purchase or our customer service (P.10).

◆ How to disposal

- When you dispose this product, please dispose it in accordance with your local regulations for sorting.
- This product uses lithium ion rechargeable battery.

Q&A

| | Question | Answer |
|-----|---|--|
| Q1 | Does it affect my eyes? | No. It does not directly affect your eyes. Please do not look directly into LED light as it may feel dark temporarily when you see a very strong light momentarily. |
| Q2 | Can I use it if I have acne? | Do not use if the inflammation is severe. When the inflammation has subsided, please restart while observing the situation gradually. |
| Q3 | Can I use it if I am tingling right after sunburn? | Please do not use it for more than three days when your epidermis is burned by ultraviolet light. You can use it if there is no stiffness or any other trouble with your skin after cleansing your face. |
| Q4 | Can I use it in other countries? | Yes. This product is compatible with 100-240V, please check the voltage of other countries. If you can change the shape of an outlet, you do not need to have a voltage converter. |
| Q5 | How long should I use this ? | It is about 24 minutes a day as a reference. It is recommended that you continue the treatment every day rather than use for a long time at a time. |
| Q6 | Can I use this in a bath? | This product is not waterproof, so please don't use it in bath, shower, or wash it with water. *Please do not store this product with water or in a humid place. |
| Q7 | Can I use this with makeup on? | Use this product after removing makeup (after cleansing and washing your face). Use this product while moisturizing your skin with lotion or beauty essences, etc. |
| Q8 | I felt face or hand is tingling. | It is electrical stimulation by EMS/desincrustation/RF/iontophoresis. There is nothing wrong with this product. If you feel the stimulation is strong, please set the function in weak level. |
| Q9 | Can I purchase spare parts if needed? Can the device be repaired? | You can do both. Please contact the shop where you purchased the device or our Exidéal customer service. |
| Q10 | How long is the life of rechargeable battery? | The battery life is about 3 years. (The warranty period of the main body is one year from the date of purchase, but consumable items including rechargeable battery are not covered by the warranty.) |

After-purchase Service

■ Warranty

The warranty is attached to the back cover. Please fill in all necessary information such as location purchased and the date of purchase and keep it safe. Please read this warranty description carefully.

■ Warranty period

The guaranteed term is one year from the date of purchase.

■ Request for repair

Please contact our customer service if you still have problems after troubleshooting.

※Shipping costs shall be borne by the customer. We pay shipping costs after repair.

● During warranty period

We will repair without charge according to the regulations of warranty. Please bring or send the product with this warranty and the proof-of-purchase, which can verify the date of purchase.

● After warranty period

When the product still performs after repairing, we will repair with charge based on the needs of customers. Please contact the original distributor or our customer service.

■ If any unsure about after service

Please contact our customer service.

Specifications

MADE IN JAPAN

| | | | |
|-----------------------|--|------------------------------------|--------------------------------|
| Battery specification | Lithium ion rechargeable battery 3.7 v 1100mAh | Approximate charging time | About 120 min |
| Power Dissipation | 2.5W | LED Output Wavelength | 468 nm, 591 nm, 632 nm, 940 nm |
| Weight | About 210 g | Number of charges | About 300 times |
| Product Dimension | Main body: about W85×H55×D121mm Charging stand: about W100×H80×D97mm Cotton ring: About W83×H7×D120mm Cable length: About1.5m | Operating temperature and humidity | 10~35°C、 10~85%RH |
| Rated power | AC100-240V 50/60Hz 0.2A | | |

* This product is a beauty device and no effects of medical equipment.

Haslux Co., Ltd.

● Beauty Business Division

ZIP 106-0032
29F ARK HILLS SENGOKUYAMA MORI TOWER, 1-9-10, Roppongi, Minato-ku, Tokyo
TEL: +81 03-6206-6188 FAX: +81 03-6206-6289

● Osaka Merchandise Center

ZIP 566-0055
2-2-6 Shinzaika, Settsu-city, Osaka
TEL : +81 06-6827-0307 FAX : +81 06-6827-0304