THANK YOU.

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ELECTRON  
ONLINE SHOP

# Safety Precautions

- Before using this product, please read the "Safety Precautions" carefully to ensure proper use.
- In the event of an abnormality, discontinue use immediately and request repair.
- Do not use the product in a malfunctioning state, as it may cause physical harm to the user.
- The following are classified as precautions and their contents are defined as follows. Important safety information is provided, so please be sure to observe it.

## 【Description of indications】

 <b>Danger</b>	Not following the contents of this indicator may result in serious injury or even the risk of death.
 <b>Caution</b>	Not following the contents of this indicator may cause personal injuries death or property damage.
 <b>Attention</b>	Not following the contents of this indicator may cause personal injuries or property damage.
 <b>Prohibited</b>	This indicator means "Prohibited" (do not do).
 <b>Mandatory</b>	This indicator means "Mandatory" (must do).

 <b>Danger</b> (Regarding use)	
 <b>Prohibited</b>	<ul style="list-style-type: none"> <li>• Do not use together with the following types of medical electronic devices. Implantable medical electronic devices such as pacemakers, life-supporting medical electronic devices such as artificial heart lungs, and wearable medical electronic devices such as electrocardiographs.</li> </ul>

 <b>Caution</b> (Regarding use)	
 <b>Prohibited</b>	<ul style="list-style-type: none"> <li>• Do not use this product or the AC adapter in a bathroom/sauna. This may cause electric shock, fire, or short circuiting.</li> <li>• Do not use this product if you fall under any of the following conditions. This may cause accidents or problems. If you have any questions, please consult your physician and follow their instructions. Acute disease/internal disease/fever/abnormal blood pressure/ infectious disease/tuberculosis/malignant tumor/blood disease/skin disease/atopic dermatitis/allergic diathesis/particularly sensitive skin/those with capillary vasodilation due to long-term use of steroidal hormones or liver dysfunction/severe fatigue/those who are not in good physical condition/those drinking alcohol or taking medications/those with sunburn irritation/those who have been prohibited by a physician, are pregnant, or menstruating.</li> </ul>

## Caution (Regarding use)

Prohibited

- Do not use on the following areas. It may result in accidents or skin/body problems.  
Areas where orthopedic surgery has been performed (It can be used if there are no problems after consulting a doctor.), inflamed areas such as eczema, rashes, acne, inflammation caused by physical irritation, etc. inside the mouth, on the eyes, throat, eyelids or ears, areas with implanted metal, plastic, silicon, etc.
- Do not look directly into the blue and red rays emitted by this product. It may result in accidents or problems.
- If abnormal heat is generated during use or charging, discontinue use immediately. Failure to do so may cause malfunction or fire.
- Do not use for any purpose other than its intended use.

## Caution (Regarding the main unit)

Prohibited

- Do not charge or store this product next to fire, close to heat, or at high temperatures. It may cause overheating, ignition, or rupture.
- Do not disassemble or modify the product. This may cause fire, electric shock, or injury.
- Store in a safe place out of reach of children. Failure to do so may result in accidents or other problems.
- Do not drop or hit this product. The main unit may malfunction, or accidents or problems may occur.
- Do not use the product with any voltage other than 100-240V. This may cause fire, electric shock, or malfunction.
- Do not damage, forcibly bend, pull, twist, bundle, place heavy objects on, or pinch the AC adapter and charging cable. This may cause fire, electric shock, or malfunction.
- Do not use the AC adapter or charging cable if it is damaged. This may cause electric shock, short circuiting, or fire.
- Do not plug or unplug the AC adapter from the power outlet with wet hands. This may cause electric shock, short circuiting, or fire.

## Caution (Regarding use)

**Mandatory**  
(Must be followed)

- If you experience pain or discomfort, discontinue use immediately.
- Use with caution to avoid skin or scalp irritation. If you experience any abnormalities or if the product does not suit you, please discontinue use. Failure to do so may cause problems.
- Those who are just starting to use the product or have sensitive skin should shorten the time of use while monitoring the condition of your skin.
- Please be sure to follow the instructions in the instruction manual for the duration and frequency of use.

## Caution (Regarding the main unit)

**Mandatory**  
(Must be followed)

- When cleaning after use, be sure to turn off the power and unplug the AC adapter from the outlet. Failure to do so may cause fire or injury.
- If an abnormality occurs during use, turn off the power, unplug the AC adapter from the outlet, and discontinue use immediately.
- Afterwards, please contact the [ELECTRON Customer Support Center] to request repair or inspection.

## Attention (Regarding use)

Prohibited

- Do not use oil or substances containing oil. This can significantly accelerate the deterioration of the comb, etc.
- Do not pull or apply strong pressure to the comb portion. This may result in damage to the product, such as the comb breaking.
- This product is not waterproof. Do not use the product in such a way that water can get inside the body except for the comb part.

## Attention (Regarding the main unit)

Prohibited

- Do not store in a humid or dusty place. This may cause malfunction.
- When cleaning, do not wipe the product with benzene/acetone/thinner. This may cause malfunction or alter the surface of the product.

## Attention (Regarding the main unit)

**Mandatory**  
(Must be followed)

- Please turn off the power to the main unit when charging. Failure to do so may cause malfunction.
- Please be sure to use the accessories specified by us for this product. Failure to do so may cause malfunction.
- When unplugging the AC adapter from the power outlet, be sure to hold the AC adapter by its sides. Failure to do so may cause electric shock, short circuiting, or fire.

## Included items

- Main unit
- Face attachment
- Body attachment
- Dedicated charging AC adapter
- Dedicated charging cable (USB Type-A to Type-C)
- Instruction manual (with warranty)
- Lotion for DENKI BARI BRUSH

### LOTION FOR DENKI BARI BRUSH

Lotion for DENKI BARI BRUSH (lotion for head, face, and entire body)

Lotion exclusively for the DENKI BARI BRUSH®.

Can be used on the entire body, including the scalp and face area.


Content: 200ml

Ingredients: Water, fullerene, rosemary leaf extract, Houttuynia cordata, Artemisia princeps extract, loquat leaf extract, repagermanium, mineral salts, BG, PVP


## ● Part names

### Face attachment

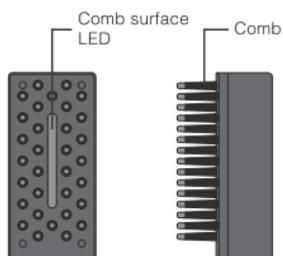It is exclusively for areas above the décolleté. Use around the scalp, face, and décolleté.

### Body attachment

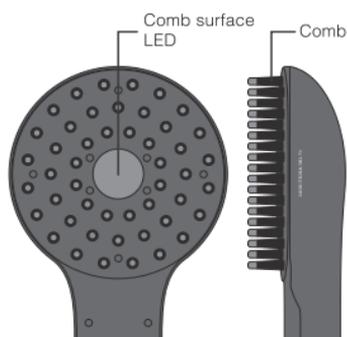
For body use only. Please use on the entire leg and abdomen area.

\* A magnet system is used to secure the attachment for easy attachment and removal.

### Main unit

## ● Main unit button/indicator names
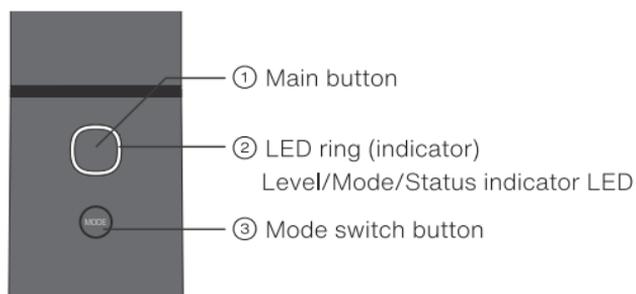
# Operation method

\* Refer to the previous section for button and indicator names.

## ● Basic operation method

- When the ① main button is pressed and held down, the ② LED ring will start blinking in white and the LED (red/blue) in the center of the comb surface will start blinking at the same time to activate the device. A notification tone (beep) will then sound indicating level 1.
- For safety reasons, the power cannot be turned on when the attachment is not attached.
- Each press of the ① main button switches the stimulation level in five steps.
- Pressing the ③ mode switch button while the main unit power is on switches the mode.
- After five minutes of continuous operation, three notification tones (beep-beep-beep) will sound, and after ten minutes, it will automatically stop with five notification tones (beep-beep-beep).
- To turn off the power during operation, press and hold the ① main button.
- This product is battery operated. When the battery power is low, the ② LED ring will flash red and a notification tone (beep) will sound continuously. Charge the battery as soon as possible.
- The ② LED ring flickers white while charging. After charging is complete, it will switch to a white lit state and turn off after a certain period of time.

## ● Switching modes

- Press the ③ mode switch button to switch modes. The color and blinking pattern of the ② LED ring changes according to the mode.
- The mode will vary depending on the attachment being worn.
- When an attachment is replaced, the main unit automatically recognizes whether the body or face is attached.

### 【Face attachment】

- There are three modes.
  - Ⓐ **Normal mode (white)**: Low-frequency stimulation is output continuously.
  - Ⓑ **Elasticity Treatment mode (purple)**: Low-frequency stimulation is output intermittently.
  - Ⓒ **Wave mode (blue)**: Low-frequency stimulation is output in waves.

### 【Body attachment】

- There are two modes.
  - Ⓐ **Normal mode (white)**: Low-frequency stimulation is output continuously.
  - Ⓑ **Elasticity Treatment mode (purple)**: Low-frequency stimulation is output intermittently.


## ● Level adjustment

- There are five levels, starting with level 1 when the power is turned on.
- The levels will switch from 1→2→3→4→5→1 with each press of the ① main button.
- The blinking of the LED ring ② will become faster as the levels increase.
- Start using with level 1.

\*Switching modes will return the level to 1.

## ● How to replace attachments


- Attachments should be replaced while the power of the main unit is turned off.
- When removing an attachment, grasp the main unit firmly and pull the attachment upward.
- When attaching an attachment, be careful not to get your hand or fingers caught in it and check the direction of the attachment carefully before inserting it straight back into the main unit.
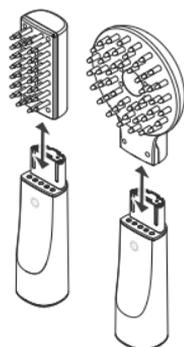
# How to use

## When using the "Face attachment"

### Scalp

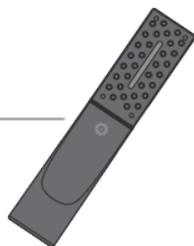- Use for about 10-20 minutes at a time, up to a maximum of 20 minutes / 2 times per day.
- Use on dry and clean scalp.
- The stimulation may vary from person to person. Please use at a level that feels comfortable to you.
- Use the Lotion for DENKI BARI BRUSH each time if you feel a strong irritation during use.
- \* An operating noise will be emitted from the product during use, but this is not a malfunction.

### Recommended modes: Normal (white)/Wave (blue)

- ① Spray the Lotion for DENKI BARI BRUSH all over the scalp and blend it in.

- ② Apply it starting from the top of the head.


- ③ Apply it along the frontal region of head.

- ④ Apply it along the area at top of head.

- ⑤ Apply it along the back of the head.

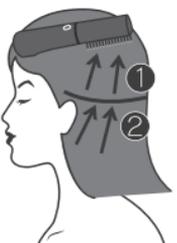
- ⑥ Move it so that it flows from the neckline to the collarbone.

- \* Since the skin around the neck is more sensitive to irritation than the scalp, please use your hands to blend the Lotion for DENKI BARI BRUSH in after spraying it on.

# How to use

## When using the "Face attachment"

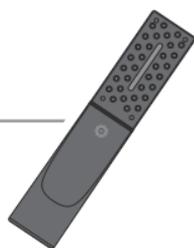
### Around the face

- Use for about 5 minutes at a time, up to a maximum of 5 minutes x 2 times per day.
  - Use on clean skin.
  - Be sure to moisturize using the Lotion for DENKI BARI BRUSH.
  - Use the Lotion for DENKI BARI BRUSH each time if you feel a strong irritation during use.
- \* An operating noise will be emitted from the product during use, but this is not a malfunction.

### Recommended modes: Normal (white)/Elasticity Treatment (purple)

- ① Spray and blend the Lotion for DENKI BARI BRUSH into the face along the neckline, neck, collarbone and face line.

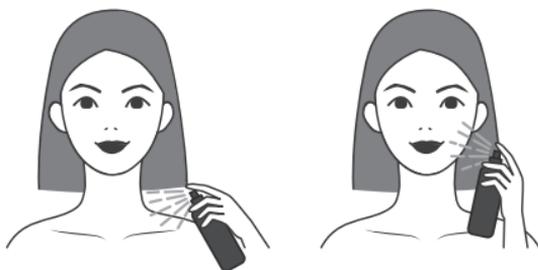- ② Move it so that it flows from the neckline to the collarbone.

- ③ Move it from under the chin, through the back of the face line, and behind the ears.

- ④ Move along the entire face in a downward and upward motion.

- ⑤ Stop moving it when you want to concentrate on a particular area of the body.

- ⑥ Finally, repeat step (2) again.

\* Do not use on the eyelids or on the Adam's apple.

# How to use

## How to use the "Body attachment"

### Legs

- Use for approximately 10-20 minutes per region.
- Be sure to moisturize using the "Lotion for DENKI BARI BRUSH" .
- Use the Lotion for DENKI BARI BRUSH each time if you feel a strong irritation during use.
- \* An operating noise will be emitted from the product during use, but this is not a malfunction.

① Spray the entire leg with the Lotion for DENKI BARI BRUSH and let it blend in.


② Move it from the base of the toes to the ankle in a flowing motion.

③ Move it along the entire calf from the ankle to the knee in a flowing motion.

④ Move it along the entire thigh in a flowing motion from the back of the knee and around the knee to the root of the leg.

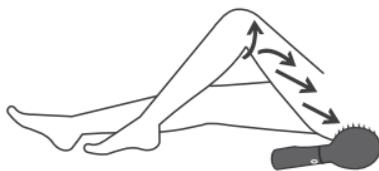SURPRISE & HAPPINESS FOR EVERYONE.

# How to use

## How to use the “Body attachment”

### Abdomen

- Use for approximately 10-20 minutes per region.
  - Be sure to moisturize using the “Lotion for DENKI BARI BRUSH” .
  - Use the Lotion for DENKI BARI BRUSH each time if you feel a strong irritation during use.
- \* An operating noise will be emitted from the product during use, but this is not a malfunction.


① Spray the Lotion for DENKI BARI BRUSH over the entire area from the tummy to the waist area and let it blend in.


② Move it in a clockwise direction around the navel.

③ Move it up and down around the waist and underarms.

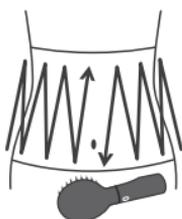④ Move it along the pelvis in a flowing motion toward the base of the foot.

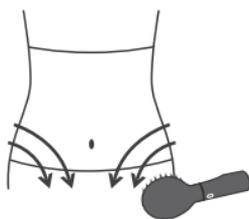It can also be used on the back, arms, hips, and other parts of the body.

\*For details, please check the “how to use” video.

A “how to use” video

is now available on our official YouTube channel.

# How to store and clean after use

- ① Please be sure to turn the power of the main unit off after use.
- ② Storage location
  - Avoid high temperature and high humidity, and store in a place with low humidity.
  - Do not store in locations where the equipment may be adversely affected by atmospheric pressure, sunlight, dusty air, etc.
  - Pay attention to safety conditions such as tilting, vibration, shock, etc.
  - Do not store in places where chemicals are stored or gases are generated.
  - Keep the product and AC adapter clean.
  - When not in use for a long period of time, it is recommended to keep the battery fully charged.
- ③ How to clean after use
  - Lightly wipe the surfaces of the equipment with a dry cloth to prevent moisture and liquids from getting inside. It cannot be sterilized with a sterilizer (e.g., UV disinfectant).
  - **After each use, gently wipe each silicon comb from base to tip with a cloth dampened with rubbing ethanol and remove any oil or hair care products before storing. Adhesion of silicone and oil components may lead to the deterioration of the comb.**
  - \* The silicon combs are consumable items. Please contact the ELECTRON Customer Support Center if you wish to replace the comb.

## ● Specifications

Power supply	AC100-240V 50/60Hz DC5V 2.4A
Battery specifications	Lithium ion battery (3.6V 960mAh 3.456Wh)
External dimensions (Including the comb portion)	Approx. W42 x H200 x D42mm (with face attachment) Approx. W100 x H241 x D42mm (with body attachment)
Weight	With face attachment: Approx. 197g With face attachment: Approx. 278g
Materials	Main unit: ABS resin Comb: Silicon (Tip: Stainless steel)
Country of manufacture	JAPAN

- Estimated continuous use time is approximately 300 minutes.  
(Depending on battery condition)
- When disposing of the product, please follow the disposal method of your area.

## ● FAQs & Troubleshooting

- If you have any problems or feel something is wrong during use, please refer to the FAQs and troubleshooting section of the official online store. If the problems persist after checking, please contact the ELECTRON Customer Support Center listed on the next page.
- If you would like to request maintenance, please apply for after-sales support for the DENKI BARI BRUSH using the QR code on the next page.

ELECTRON  
ONLINE SHOP

# Warranty provisions

---

In the event of a malfunction during the following period from the date of purchase, this document warrants to repair the product free of charge in accordance with the contents described in this document.

The use of the personal information provided by you is described below.

1. If the product malfunctions within the warranty period under conditions of use in accordance with the precautions set forth in the instruction manual, etc., the main unit and attachments (excluding the comb portion) will be repaired free of charge for one year from the date of purchase, and the comb, dedicated charging AC adapter and dedicated charging cable will be repaired free of charge for six months.
  2. Repairs will be charged in the following cases even if it is within the warranty period.
    - (1) Malfunction or damage due to incorrect use, improper repair or modification.
    - (2) Malfunction or damage due to transportation, dropping, etc. after purchase
    - (3) Malfunction or damage due to fire, earthquake, flood, lightning or other natural disasters, pollution, salt damage, gas damage, or abnormal voltage
    - (4) Malfunction or damage caused by assembling on a vehicle, ship, etc.
    - (5) When replacing batteries or worn parts.
    - (6) When purchased as a second-hand product, through private sales, online auctions, or resale stores on the Internet.
  3. This document is only valid in Japan.
  4. This document will not be reissued, so please keep it in a safe place along with an itemized invoice showing the date of purchase.
- \* Please note that your personal information (copy of the warranty) may be used to provide free repair services during the warranty period and for subsequent safety inspection activities.

**Manufacturer: GM CORPORATION**

3-6-8 Toyosaki, Kita-ku, Osaka-shi, Osaka Prefecture, 531-0072

**ELECTRON Customer Support Center**

Toll-free number

**0120-66-7170**

Reception hours: Monday-Friday 11:00-18:00 (excluding national holidays)

**DENKI BARI BRUSH After-sales Service**

<https://gm-beauty.jp/contact/>

\* Scan the QR code to activate the input form.

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