THANK YOU.

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Safety Precautions

- Before using this product, please read the "Safety Precautions" carefully to ensure proper use.
- In the event of an abnormality, discontinue use immediately and request repair.
- Do not use the product in a malfunctioning state, as it may cause physical harm to the user.
- The following are classified as precautions and their contents are defined as follows. Important safety information is provided, so please be sure to observe it.

【Description of indications】

 Danger	Not following the contents of this indicator may result in serious injury or even the risk of death.
 Caution	Not following the contents of this indicator may cause personal injuries death or property damage.
 Attention	Not following the contents of this indicator may cause personal injuries or property damage.
 Prohibited	This indicator means "Prohibited" (do not do).
 Mandatory	This indicator means "Mandatory" (must do).

 Danger (Regarding use)	
 Prohibited	• Do not use together with the following types of medical electronic devices. Implantable medical electronic devices such as pacemakers, life-supporting medical electronic devices such as artificial heart lungs, and wearable medical electronic devices such as electrocardiographs.

Caution (Regarding use)

Prohibited

- Do not use this product or the AC adapter in a bathroom/sauna. This may cause electric shock, fire, or short circuiting.
- Do not use this product if you fall under any of the following conditions. This may cause accidents or problems. If you have any questions, please consult your physician and follow their instructions. Acute disease/internal disease/fever/abnormal blood pressure/infectious disease/tuberculosis/malignant tumor/blood disease/skin disease/atopic dermatitis/allergic diathesis/particularly sensitive skin/those with capillary vasodilation due to long-term use of steroidal hormones or liver dysfunction/severe fatigue/those who are not in good physical condition/those drinking alcohol or taking medications/those with sunburn irritation/those who have been prohibited by a physician, are pregnant, or menstruating.
- Do not use on the following areas. It may result in accidents or skin/body problems.
Areas where orthopedic surgery has been performed, inflamed areas such as mucous membranes, eczema, throat sores, rashes, acne, etc. inside the mouth, on the eyes, eyelids or ears, areas with implanted metal, plastic, silicon, etc., itching, inflammation caused by physical irritation, stains (pigmentation), etc.
- Do not look directly into the blue and red rays emitted by this product. It may result in accidents or problems.
- If abnormal heat is generated during use or charging, discontinue use immediately. Failure to do so may cause malfunction or fire.
- Do not use for any purpose other than its intended use.

Caution (Regarding the main unit)

Prohibited

- Do not charge or store this product next to fire, close to heat, or at high temperatures. It may cause overheating, ignition, or rupture.
- Do not disassemble or modify the product. This may cause fire, electric shock, or injury.
- Store in a safe place out of reach of children. Failure to do so may result in accidents or other problems.

Caution (Regarding the main unit)

Prohibited

- Do not drop or hit this product. The main unit may malfunction, or accidents or problems may occur.
- Do not use the product with any voltage other than 100-240V. This may cause fire, electric shock, or malfunction.
- Do not damage, forcibly bend, pull, twist, bundle, place heavy objects on, or pinch the AC adapter and charging cable. This may cause fire, electric shock, or malfunction.
- Do not use the AC adapter or charging cable if it is damaged. This may cause electric shock, short circuiting, or fire.
- Do not plug or unplug the AC adapter from the power outlet with wet hands. This may cause electric shock, short circuiting, or fire.

Caution (Regarding use)

Mandatory (Must be followed)

- If you experience pain or discomfort, discontinue use immediately.
- Use with caution to avoid skin or scalp irritation. If you experience any abnormalities or if the product does not suit you, please discontinue use. Failure to do so may cause problems.
- Those who are just starting to use the product or have sensitive skin should shorten the time of use while monitoring the condition of your skin.
- Please be sure to follow the instructions in the instruction manual for the duration and frequency of use.

Caution (Regarding the main unit)

Mandatory (Must be followed)

- When cleaning after use, be sure to turn off the power and unplug the AC adapter from the outlet. Failure to do so may cause fire or injury.
- If an abnormality occurs during use, turn off the power, unplug the AC adapter from the outlet, and discontinue use immediately.
- Afterwards, please contact the [ELECTRON Customer Support Center] to request repair or inspection.

Attention (Regarding use)

Prohibited

- Do not use oil or substances containing oil. This can significantly accelerate the deterioration of the comb, etc.
- Do not pull or apply strong pressure to the comb portion. This may result in damage to the product, such as the comb breaking.
- This attachment is for body use only. Do not use it around the head & face. Failure to do so may cause problems.

Attention (Regarding the main unit)

Prohibited

- Do not store in a humid or dusty place. This may cause malfunction.
- When cleaning, do not wipe the product with benzene/acetone/thinner. This may cause malfunction or alter the surface of the product.

Attention (Regarding the main unit)

Mandatory
(Must be followed)


- Please turn off the power to the main unit when charging. Failure to do so may cause malfunction.
- Please be sure to use the accessories specified by us for this product. Failure to do so may cause malfunction.
- When unplugging the AC adapter from the power outlet, be sure to hold the AC adapter by its sides. Failure to do so may cause electric shock, short circuiting, or fire.

SURPRISE & HAPPINESS FOR EVERYONE.


Included items

Body attachment

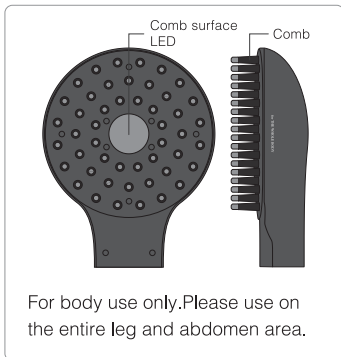
Instruction manual (with warranty)

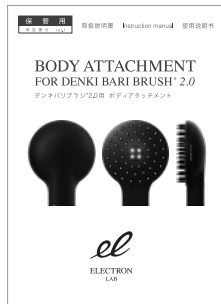*A magnet system is used to secure the attachment for easy attachment and removal.

● How to replace attachments

- Attachments should be replaced while the power of the main unit is turned off.
- When removing an attachment, grasp the main unit firmly and pull the attachment upward.
- When attaching an attachment, be careful not to get your hand or fingers caught in it and check the direction of the attachment carefully before inserting it straight back into the main unit.

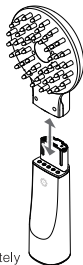*Main unit: Sold separately

● Specifications

Dimensions	Approx. W100×H127×D42mm
Weight	137g
Material	Main unit: ABS resin Comb: Silicon (Tip: Stainless steel)
Country of manufacture	JAPAN

*When disposing of the product, please follow the disposal method of your area.

Warranty provisions

In the event of a malfunction during the following period from the date of purchase, this document warrants to repair the product free of charge in accordance with the contents described in this document.

The use of the personal information provided by you is described below.

1. If the product malfunctions within the warranty period under conditions of use in accordance with the precautions set forth in the instruction manual, etc., the main unit and attachments (excluding the comb portion) will be repaired free of charge for one year from the date of purchase, and the comb, dedicated charging AC adapter and dedicated charging cable will be repaired free of charge for six months.
2. Repairs will be charged in the following cases even if it is within the warranty period.
 - (1) Malfunction or damage due to incorrect use, improper repair or modification.
 - (2) Malfunction or damage due to transportation, dropping, etc. after purchase
 - (3) Malfunction or damage due to fire, earthquake, flood, lightning or other natural disasters, pollution, salt damage, gas damage, or abnormal voltage
 - (4) Malfunction or damage caused by assembling on a vehicle, ship, etc.
 - (5) When replacing batteries or worn parts.
 - (6) When purchased as a second-hand product, through private sales, online auctions, or resale stores on the Internet.
3. This document is only valid in Japan.
4. This document will not be reissued, so please keep it in a safe place along with an itemized invoice showing the date of purchase.

* Please note that your personal information (copy of the warranty) may be used to provide free repair services during the warranty period and for subsequent safety inspection activities.

Manufacturer: GM CORPORATION

3-6-8 Toyosaki, Kita-ku, Osaka-shi, Osaka Prefecture, 531-0072

 **ELECTRON Customer Support Center**

Toll-free number **0120-66-7170**

Reception hours:

Monday-Friday 11:00-18:00 (excluding national holidays)

 **DENKI BARI BRUSH After-sales Service**

<https://gm-beauty.jp/contact/>

* Scan the QR code
to activate the input form.

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